

Specifications:


Dimensions (L x W x H):	19 x 13 x 16in.(490 x 330 x 406mm)
Weight:	99lbs.(45 kg)
Power requirements:	AC,90-240V, 50/60 Hz, 350 Watts
Laser model:	Ultraviolet laser
Laser power:	1 Watt
Working temperature:	59°F ~ 86°F (15°C ~ 30°C)
Work noise:	≤65dB
Connectivity:	USB2.0
Cassette printing speed:	< 3 seconds per cassette (on-demand printing)
High resolution:	500 DPI
Number of hoppers :	4 hoppers for one piece tissue cassette
Cassette capacity per hopper:	80 cassettes per hoppers
Total cassette capacity:	320 cassettes
Output tray capacity:	7 cassettes

Guangdong Jinqun Medical tissue cassette laser printer is a must-have for pathology and lab environments. Precision construction forges jam-free printing while its low work noise ≤65dB design offers ribbon or cartridge-free printing, optimal for cost-efficient workflows.

Easily integrated into laboratory and hospital information systems, this laser cassette printer is perfectly compatible with our hinged tissue cassettes and able to print on cassettes with or without lids. Its structure possesses 4 cassette magazines holding up to 320 cassettes total. Choose from many printing options, including 1D as well as 2D barcoding, alphanumeric, graphic, and logo printing. On-demand mode printing, its high printing speed completes a cassette every 3 seconds.

Features:

- Laser tissue cassettes printer (ribbon or cartridge free printing)
- 4 cassette magazines hold up to 320 cassettes total, load 4 different color cassettes for efficient workflow
- Precision construction for jam-free cassette printing
- 1D and 2D barcodes, alphanumeric, graphic and logo printing
- Low work noise ≤65dB, meet the hospital pathology and labs environment requirement
- Integrate to LIS/HIS System
- Perfectly compatible with our hinged tissue cassettes
- High printing speed completes a cassette every 3 seconds

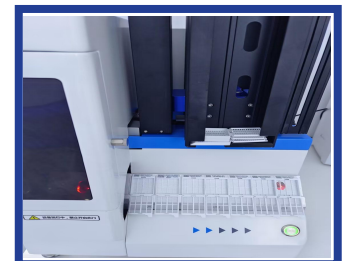

Model No.	Description
EMR840	Laser Cassette Printer